



## EMS System for Metropolitan Oklahoma City and Tulsa 2026 Medical Control Board Treatment Protocols



Approved 1/7/26, Effective 4/1/26, Replaces all prior versions

### 4E – DOUBLE SEQUENTIAL EXTERNAL DEFIBRILLATION ADULT

#### PARAMEDIC

##### Indication:

**Adult** refractory Ventricular Fibrillation/Pulseless Ventricular Tachycardia. See also Protocol 4G - Ventricular Fibrillation & Pulseless Ventricular Tachycardia.


##### Contraindications:

Spontaneous pulse.

All cardiac rhythms except ventricular fibrillation/pulseless ventricular tachycardia

##### Technique:

1. If two monitor/defibrillators are available, power the second one ON.
2. Connect therapy electrodes (defibrillation pads) to therapy cable and confirm cable connection to monitor/defibrillator per Protocol 4D.
3. Prepare the patient's skin and apply second set of therapy electrodes (defibrillation pads) to the patient, in the right parasternal and cardiac apex positions next to, but NOT overlapping the anterior pad of the first set of therapy electrodes (defibrillation pads). (Figure 1)
4. Proceed to charge each defibrillator to 360J. Once fully charged, have either one paramedic discharge both defibrillators or if using two paramedics, using a 3-2-1 verbal countdown, **discharge the defibrillators sequentially using a very specific count of "one thousand one" between the discharge of defibrillator one and the discharge of defibrillator two. FAILURE TO ALLOW A ONE SECOND PAUSE BETWEEN SEQUENTIAL DEFIBRILLATIONS CAN IRREVERSIBLY HARM A DEFIBRILLATOR, RENDERING IT PERMANENTLY INOPERABLE.**

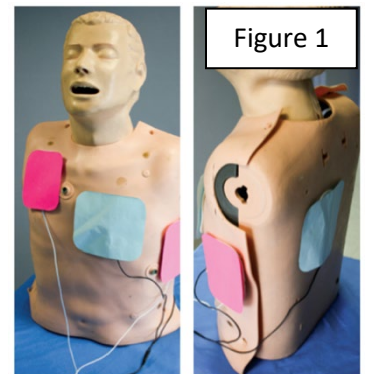



## EMS System for Metropolitan Oklahoma City and Tulsa 2026 Medical Control Board Treatment Protocols



### PROTOCOL 4E: – Double Sequential External Defibrillation - Adult, cont.

5. If only one monitor/defibrillator is available and the treatment algorithm indicates the **administration of double sequential external defibrillation**, change the defibrillation energy vector by applying a secondary set of defibrillation pads (position of secondary set of pads right parasternal and cardiac apex as outlined in step 3 of this treatment protocol). If a second monitor/defibrillator becomes available, follow the standard double sequential external defibrillation process.
6. Throughout the use of double sequential external defibrillation, follow all standard safety measures as with routine defibrillation as outlined in Protocol 4D.

#### CLINICAL PEARLS:

**DO NOT CONTINUE TO PROVIDE CHEST COMPRESSIONS WHEN THE MONITOR/DEFIBRILLATOR IS DISCHARGING / DEFIBRILLATING.**